
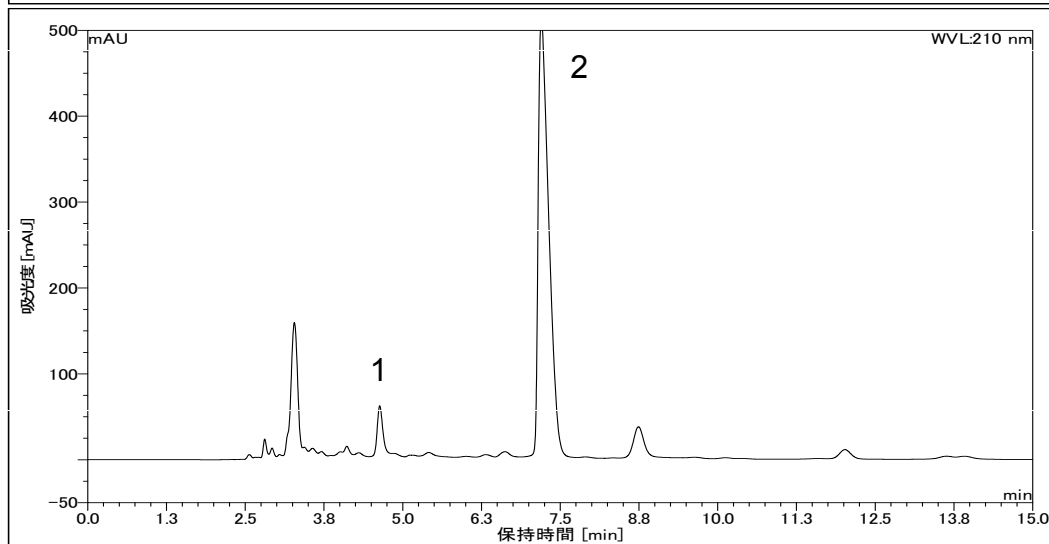


Ascorbic Acid and Citric Acid contained in Drinks

ドリンク中に含まれるアスコルビン酸・クエン酸の分析

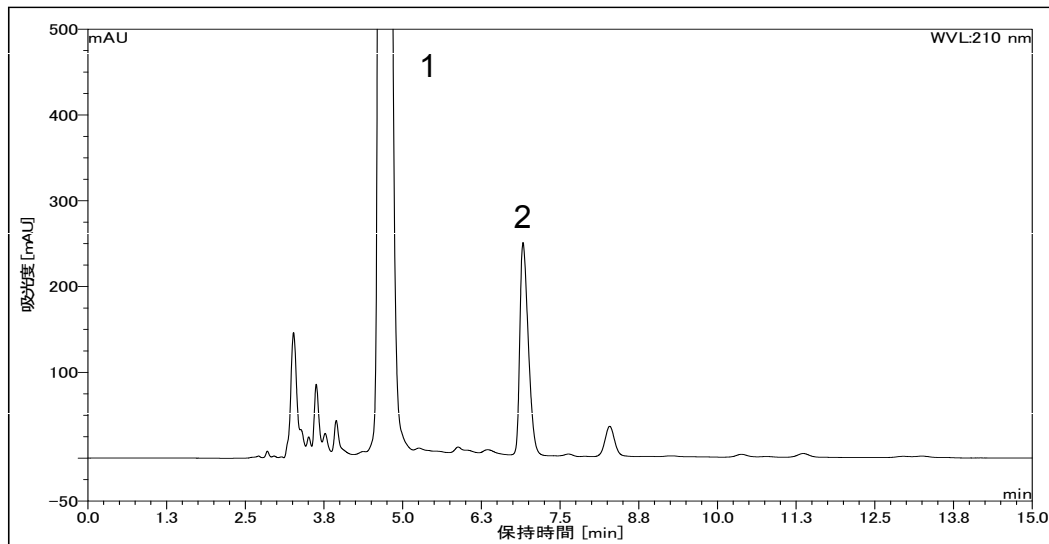

**Conditions;**

**Column:** Develosil XG-C1 (5um)

**Size:** 4.6x250mm

**Mobile phase:** Acetonitrile/0.2%Phosphoric acid=2/98

**Flow rate:** 1.0ml/min


**Temperature:** 30°C

**Detection:** UV210nm (PDA)


**Sample:** 1.Ascorbic acid 2.Citric acid

**Injection volume:** 10uL

**System:** DIONEX/Ultimate 3000 Series

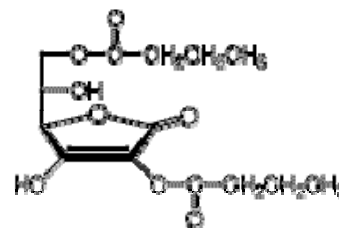

Ascorbic acid

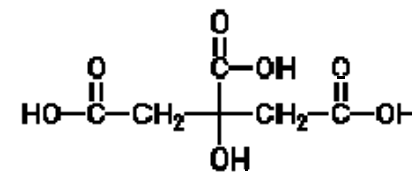

Citric acid