


Analysis of Pesticides (4)/ 農薬の分析

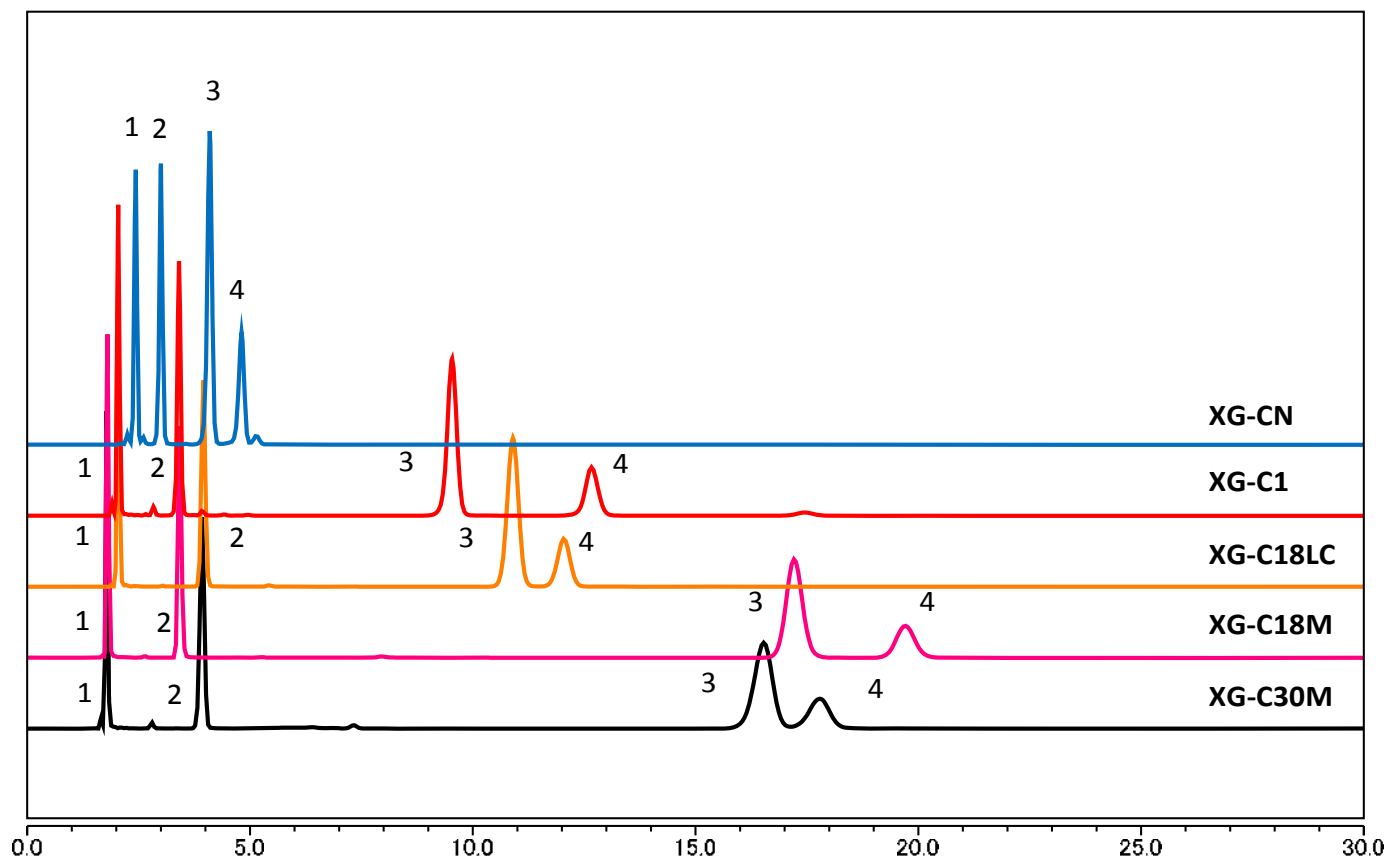

Analytical conditions	
Column	Develosil XGシリーズ(5um) 4.6x250mm
Mobile phase	Methanol/0.2%Phosphoric acid=60/40
Flow rate	1.0ml/min
Temperature	40°C
Detection	UV210nm
Sample	1: Asulam (0.44mg/mL) 2: Thiuram (0.62mg/mL) 3: Iprodione (0.49mg/mL) 4: Bensulide (0.78mg/mL)
Injection volume	5.0uL
System	Thermo Scientific ultimate 3000 series