


Develosil Application Data

No.086 Peptide/ペプチド

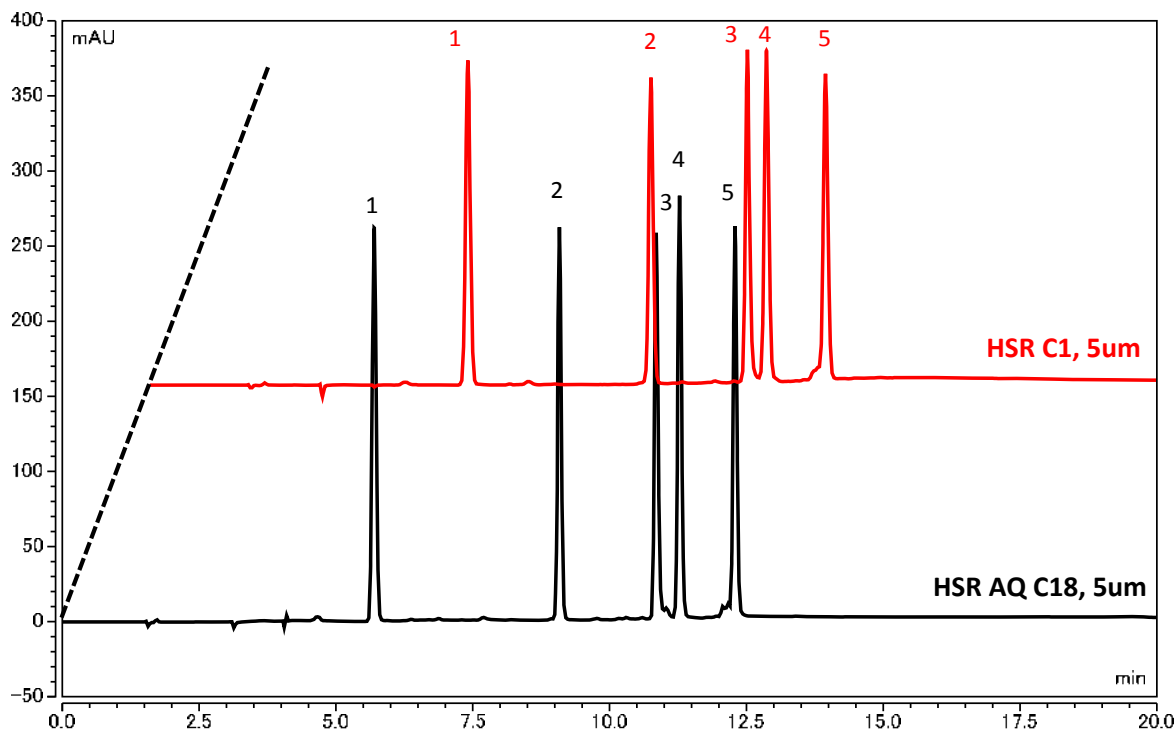

Conditions;

Column: Develosil HSR C1, 5um (4.6x150mm)

Develosil HSR AQ C18, 5um (4.6x150mm)

Mobile phase: A) 25mM NaCl (pH3.0) B)Acetonitrile

Flow rate: 1.0mL/min

Temperature: 30°C

Sample: 1.Gly-Tyr (0.05mg/mL)

2.Val-Tyr-Val (0.05mg/mL)

3.Met-Enkephalin (0.05mg/mL)

4.Angiotensin II (0.05mg/mL)

5.Leu-Enkephalin (0.05mg/mL)

Injection volume: 10.0uL

System: Thermo SCIENTIFIC Ultimate 3000

Gradient

Time (min)	Flow (mL/min)	%A	%B	Curve
0.0	1.0	100	0	5
10.0	1.0	70	30	5
20.0	1.0	70	30	5
20.1	1.0	100	0	5

