


Analysis of Nucleobases/核酸塩基の分析

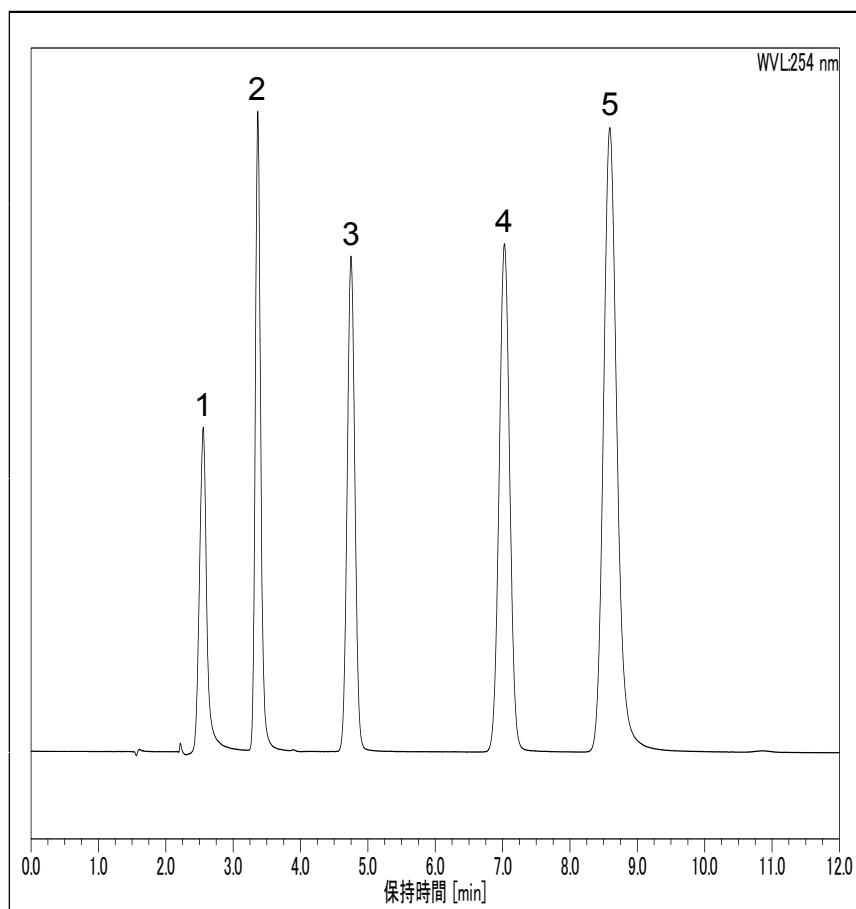

Sample

1.Cytosine (0.063mg/ml)

2.Cytidine (0.16mg/ml)

3.Uracil (0.062mg/ml)

4.Uridine (0.18mg/ml)

5.Thymine (0.20mg/ml)

Conditions;

Column: Develosil XG-C18M (5um) 4.6x150mm

Mobile phase: 50mM Phosphate Buffer (pH7.0)

Flow rate: 1.0ml/min

Temperature: 40°C

Detection: UV254nm

Injection volume: 2.0uL

System: DIONEX/Ultimate 3000 Series