



Analysis of Crospovidone/日本薬局方 クロスロビドン

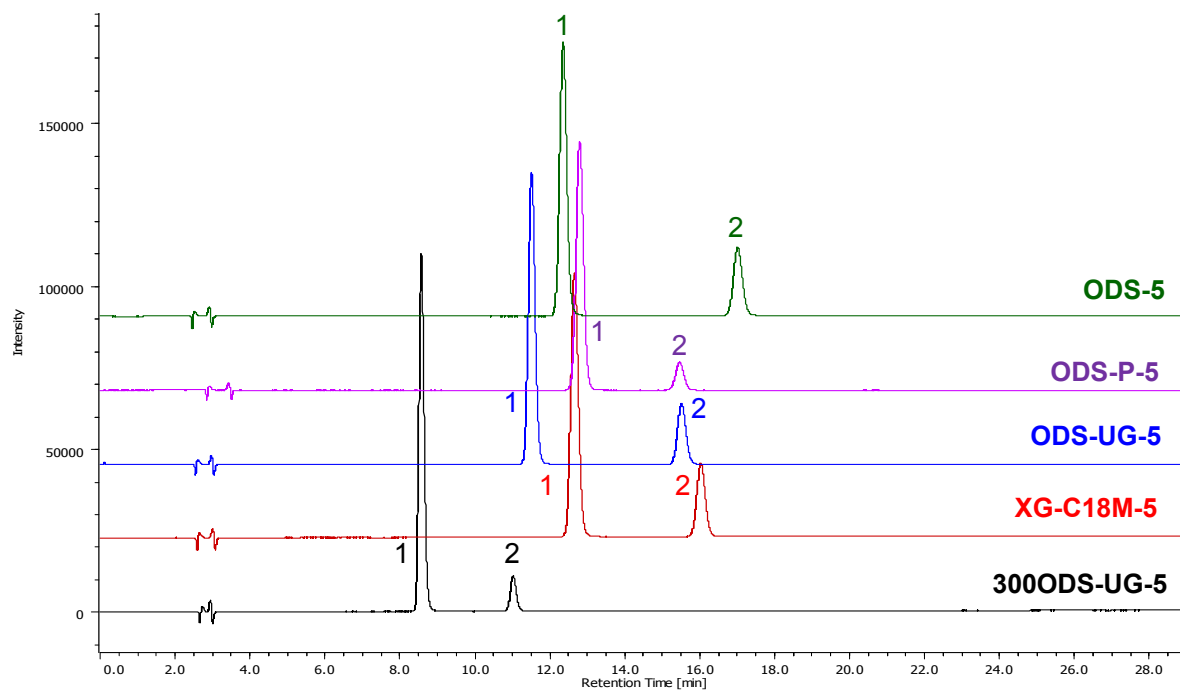


Sample:

1. 1-Vinyl-2-pyrrolidone
2. Vinyl Acetate

※サンプルの調製は局方記載に従う

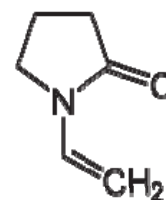

1-Vinyl-2-pyrrolidone

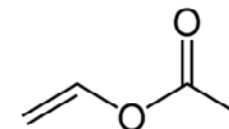

Vinyl Acetate

Conditions;

Column: Develosil 300ODS-UG-5,XG-C18M-5,ODS-UG-5,ODS-P-5,ODS-5

Size: 4.0x25mm + 4.0x250mm

Mobile phase: Acetonitrile/Water=10/90

Flow rate: 1.0ml/min

Temperature: 40°C

Detection: UV235nm

System: 日本分光/LC-2000plus series