


Develosil Application Data

No.068 Oil-soluble vitamins/脂溶性ビタミン

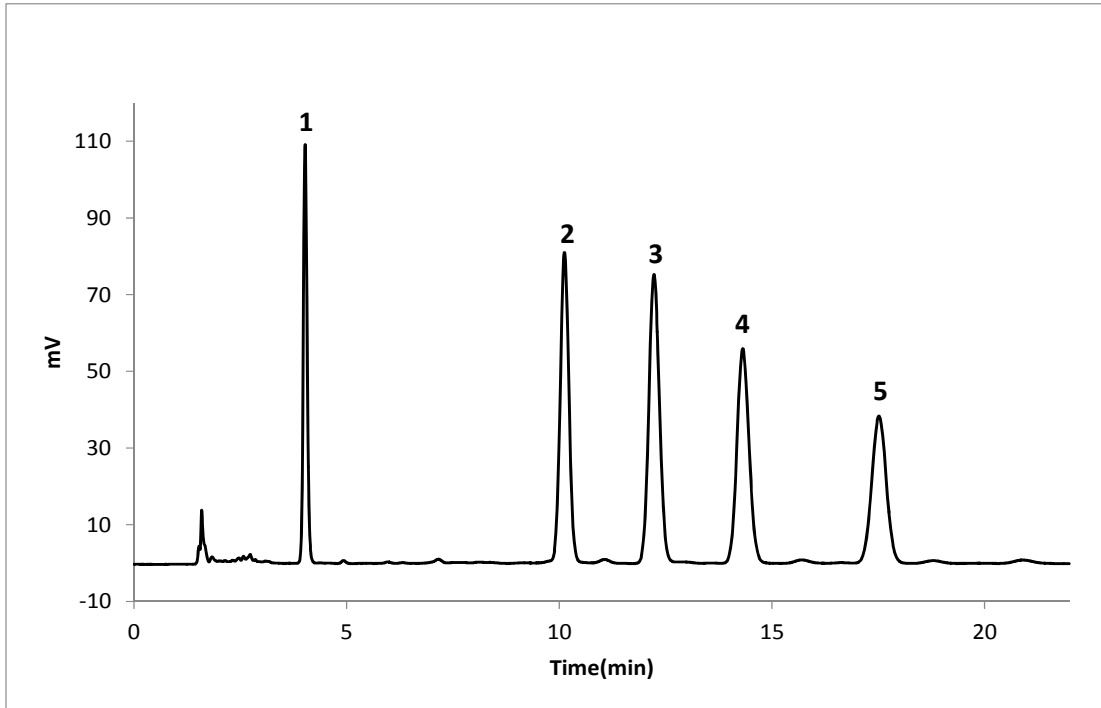

Analytical Conditions:

Column: Develosil HSR C18,5 μ m (4.6x150mm)

Mobile phase: Acetonitrile/Methanol=60/40

Flow rate: 1.0mL/min

Temperature: 30°C

Detection: UV280nm

Sample: 1.Vitamin A acetate (0.60mg/mL)

2.VitaminD3 (0.63mg/mL)

3.Vitamin E (5.21mg/mL)

4.Vitamin E acetate (5.29mg/mL)

5.Vitamin K1 (1.89mg/mL)

Injection volume: 1.0 μ L

System: JASCO2000 series

