


局方：エストリオールの分析/Analysis of Estriol

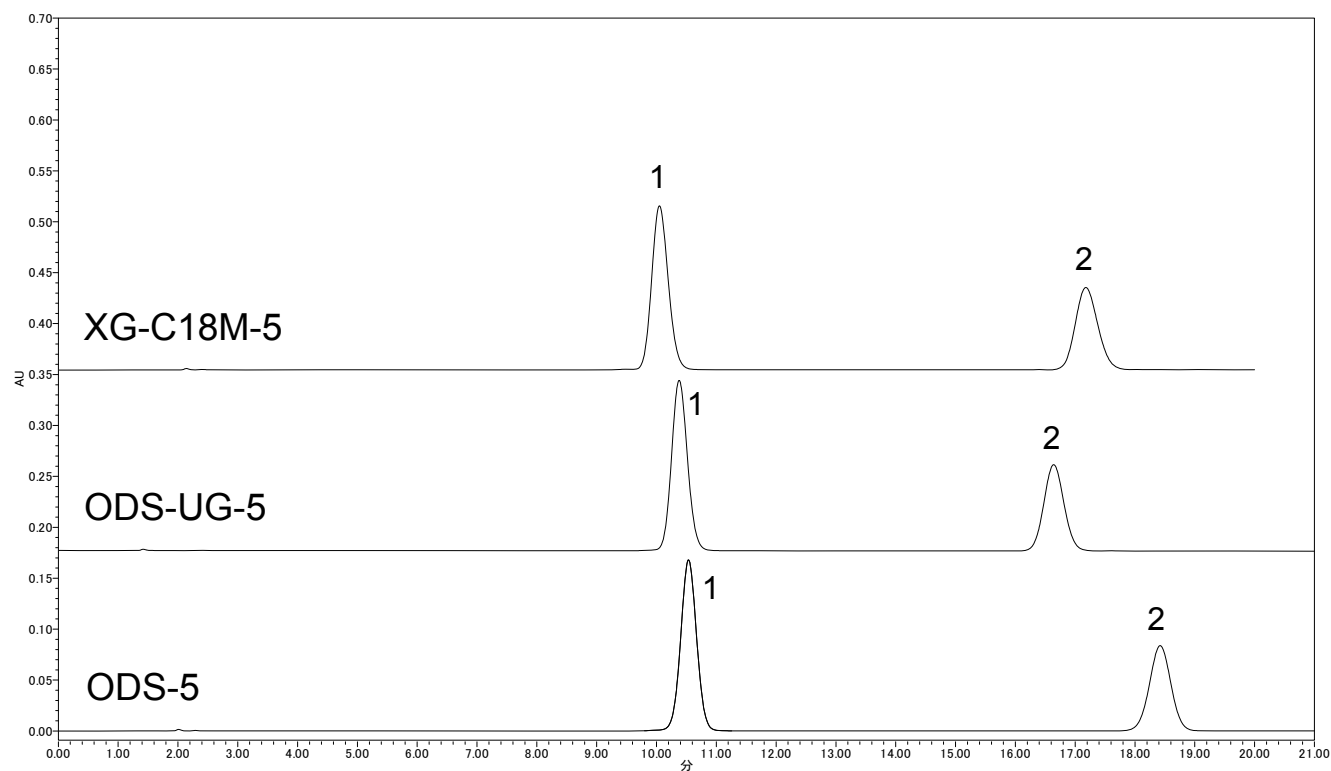

Analytical conditions	
Column	Develosil XG-C18M-5 (4.6x250mm)
	Develosil ODS-UG-5(4.6x250mm)
	Develosil ODS-5 (4.6x250mm)
Mobile phase	In accordance to Japanese pharmacopoeia for estriol
Flow rate	0.7ml/min
Temperature	27°C
Detection	UV280nm
Sample	1: Estriol 2: Methyl benzoate
Injection volume	10uL
System	Waters Alliance