



Develosil[®] Application Data No.053

Herbal Medicine Rhubarb/日本薬局方 ダイオウ


Develosil ODS-5

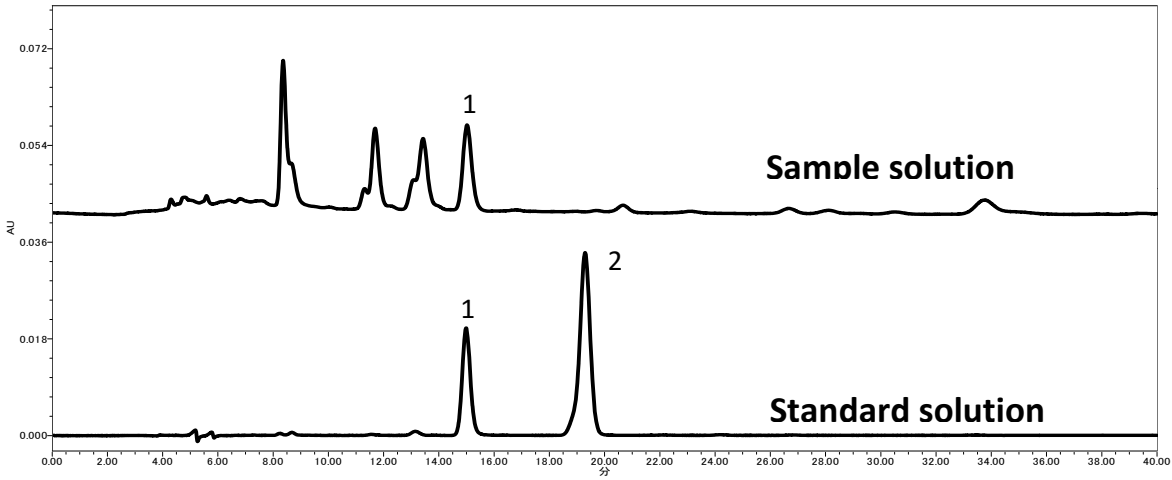

Develosil ODS-UG-5

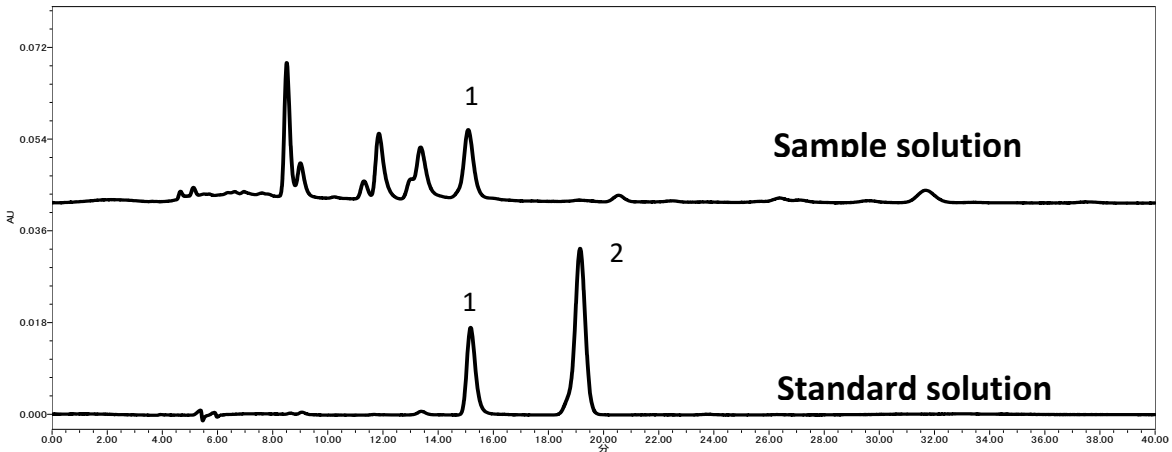

Develosil XG-C18M-5

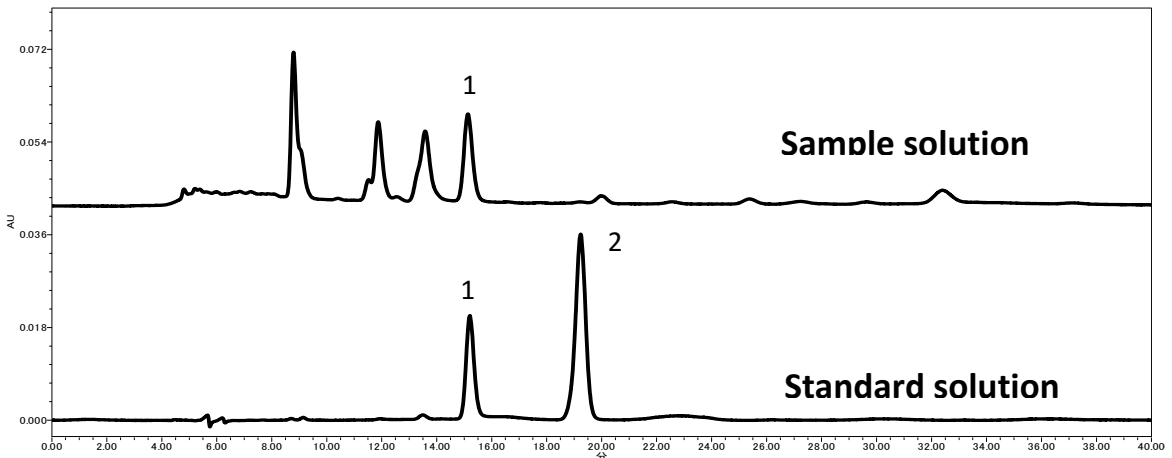

Column: Develosil ODS-5 (4.6x150mm)

Develosil ODS-UG-5 (4.6x150mm)

Develosil XG-C18M-5 (4.6x150mm)

Mobile Phase: 1.25% Acetic acid/Acetonitrile=4/1

Flow rate: 0.34mL/min (Develosil ODS-5)

0.35mL/min (Develosil ODS-UG-5)

0.33mL/min (Develosil XG-C18M-5)

Temperature: 40°C

Sample: Standard solution

1: Sennoside A (0.05mg/mL)

2: Naringin (0.05mg/mL)

Injection Volume: 20µL

Detection: UV340nm

	Requirement	Result
Resolution(1,2)	≥ 3	7.02 (Develosil ODS-5)
		6.36 (Develosil ODS-UG-5)
		6.69 (Develosil XG-C18M-5)